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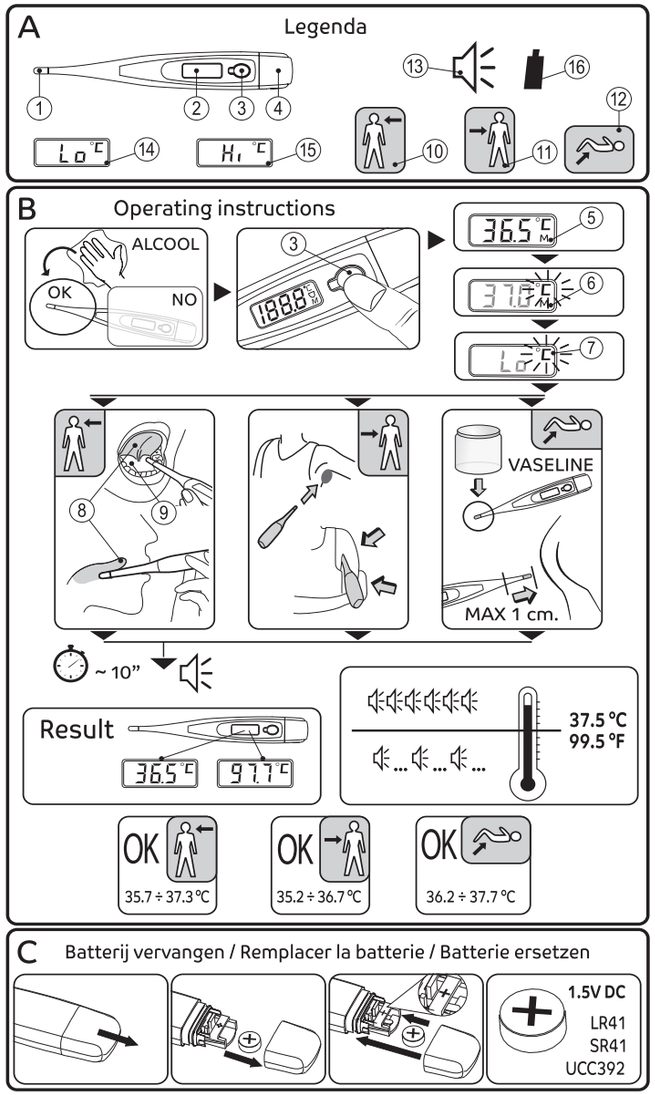
EN

Digital Thermometer

USE DESTINATION
The DIGITAL/ELECTRONIC THERMOMETER is a medical device used as an help to the medical check to keep body temperature under control. The device is destined to be used by untrained users. Incapable people, children, aged people can used it if supported by a capable adult.

LEGENDA (A)

1. Probe
2. Display
3. ON/OFF button
4. Battery Cap
5. Last measurement memory
6. Ready for measurement
7. Please wait
8. Tongue
9. Undertongue
10. Oral measurement
11. Axillary measurement
12. Rectal measurement
13. Acoustic alarm
14. Temperature is below 32 °C
15. Temperature is above 42.9 °C
16. Low battery

⚠ WARNINGS

- Failure to follow the warnings and the content of the present manual can expose to heavy risks. Read instructions carefully before using the digital thermometer.
- Never leave the appliance and all its parts, including packaging, within the reach of children, uncapable people, domestic animals or insects. The appliance includes some small parts which can be dangerous if inhaled or swallowed.
- If the thermometer has been stored at below-freezing temperatures, allow it to warm naturally to room temperature before using.
- Do not drop the thermometer or expose it to heavy shock. It is not shock-proof.
- Do not expose the thermometer to strong electromagnetic interference such as microwave ovens and cell phones.
- Do not expose to direct sunlight, high moisture or dust
- Do not dismount the device.
- Do not wet the device, clean it with a dry cloth.
- Contact always your doctor for the choice of any therapy.
- Stop using the thermometer if it operates erratically.

USE INSTRUCTIONS (B)
Disinfect the tip of the thermometer with alcool before and after use.
Perform the measurement as described in picture B.
When the “Low Battery (16)” symbol shows up, replace the battery as per pict. C.
Dispose of thermometer and batteries in compliance to regulations in force about electronic devices disposal.

WARRANTY
The thermometer is guaranteed for one year from the date of purchase (two years only for EU countries) against manufacturer’s defect under normal, household use. If the thermometer does not function properly, first check the battery.

TECHNICAL DATA
Measurement range: 32-42.9 °C
Accuracy: 0.1 °C - Resolution: ±0.1 °C
Environmental working conditions:
Temperature: 5-40 °C; Humidity: 15-95%

CE				LOT	
Compliant to european directive 93/42/EEC “0051” is the identification number of the notified body	Type BF	Attention: read carefully the documentation	Dispose of as per regulation in force	Lot number	Manufacturing date

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Manufacturer : Norditalia Group S.r.l.